


Fluoride Oxalate Tubes

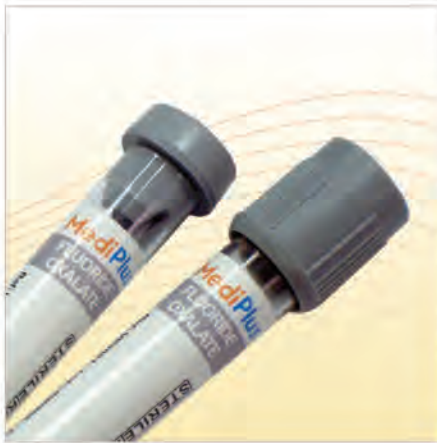

Purpose:

These tubes are generally used for glucose tests. Sodium Fluoride and Oxalate are added to decrease enzymatic activity and glucose metabolism in collected sample.

Additives:

Sodium Fluoride and Oxalate

Recommendation:

Invert and back the collected sample 8 to 10 times. For centrifugation, run the sample at $\leq 1300g$ for 10 minutes.

Specifications:

Product No.	Volume	Size(mm)	Additives	Material	Box/Case Qty
VP2091	2ml	13 x 75	Fluoride Oxalate	PET	100/1000
VP3091	3ml	13 x 75	Fluoride Oxalate	PET	100/1000

Needles & Holders



- Bevel angle is carefully calculated to provide resistance that allows easy penetration and prevents discomfort to patients.
- Advanced thin needle walls that help preserve blood cells.
- Sterilized with Ethylene Oxide.
- Holders fit perfectly with blood collection tubes and can be used repeatedly.

Product No.	Product Name	Cap Color	Specs.	Box Qty
G1810	Blood Collection Needle	Pink	18G, 1"	100
G1815	Blood Collection Needle	Pink	18G, 1.5"	100
G2015	Blood Collection Needle	Yellow	20G, 1.5"	100
G2115	Blood Collection Needle	Green	21G, 1.5"	100
G2215	Blood Collection Needle	Black	22G, 1.5"	100
H0001	Needle Holder	Clear		200



Size Chart

Our blood collection tubes can also be made into custom specifications.

Size (mm)	Volume
13 x 75	1.8ml, 2ml, 2.5ml, 2.7ml, 3ml, 3.5ml, 4ml, 4.3ml, 4.5ml, 5ml
13 x 100	4ml, 4.5ml, 5ml, 6ml, 7ml
16 x 100	8ml, 8.5ml, 9.5ml, 10ml